

NEL19-PDO-02_A05

Benefits

- Class 1 Medical Device
- Hermetic device for optimum drying
- Hermetic device for a perfect UV-C treatment
- No risk of UV exposure
- USB cable and plug included
- 2 USB connectors for simultaneous charging and drying
- User-friendly
- one size fits all chargers



- Multiple use



PerfectDome
Hear well every day

PerfectDome

Drying & UV-C System

new



Specifically designed
for chargers and hearing aids


Hear well every day

mg
Development The Hearing Hygiene Specialist

MG Développement
Zac les portes de l'aéroport
178 rue Negue-Cat
34130 Manguio - FRANCE



Easy to use

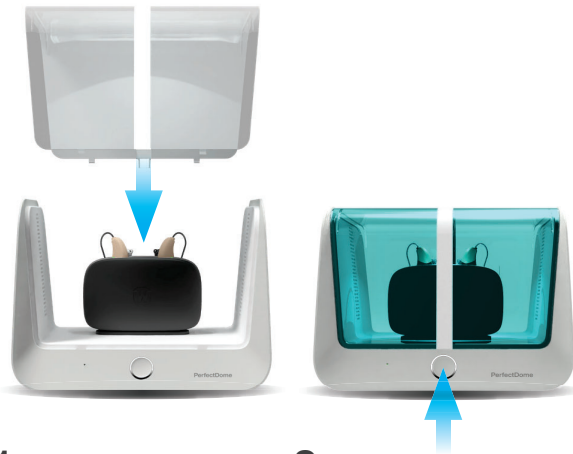

1
With its curved shape, even users lacking dexterity can easily place the lid.

2
Simply press the button to start the automatic cycle (Drying and UV-C).

Connectors

2 USB connectors (with silicone cable guide) allow to combine maintenance and charging of hearing aids, regardless of the brand.





Drying system

Humidity is the main cause of hearing aid failure.

The PerfectDome incorporates a drying system which include a fan and heater. The designed air flow combined with a sealed, hermetic chamber ensure that it dries all surfaces.

By keeping the temperature inside the device under control, drying is more efficient and has minimum impact on the battery durability. It also eliminates moisture which offer optimum hearing aid performance.

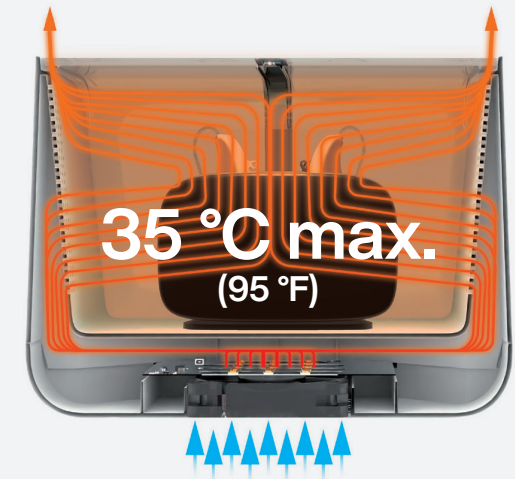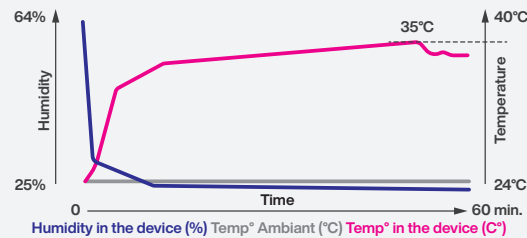

UV-C system

The use of hearing aids causes irritation and itching in the ear canal, which creates discomfort for the end-user.

The PerfectDome integrates a UV-C cleaning system based on 5 LEDs that are positioned as close as possible to the hearing aid.

A regular use of our UV-C devices will reduce significantly ear irritations, increasing customer satisfaction.

In terms of user safety, as the device is fully closed, it also prevents unwanted UV-C exposure.

